

NANOPASS™ 32.5G 4mm 6mm 8mm

Pen needle fits all common pen-injectors.⁴



Specifications

Nanopass™ 32.5G
(Ø 0.22 mm)

	4 mm	6 mm	8 mm
Needle length	4 mm	6 mm	8 mm
Product code	TN*32504M	TN*32506M	TN*32508M
Packaging size	100 pieces	100 pieces	100 pieces

CE 0197 Does not contain components made of natural rubber latex.

Visit www.terumo-europe.com

¹ Patented technology.

² Intra-company data show less puncture resistance for the asymmetric edge compared to a standard 3-bevel. Lower puncture resistance is related to lower pain (T. Asakura, 2013, Jpn J Med Pharm Sci 69(1):147-155).

³ The size of the perforation scar correlates with the outer diameter of the needle (T. Asakura, H. Seino, 2004, Jpn J Pharm. Health Care Sci., 30(6):368-376).

⁴ For a detailed list of compatible devices see NANOPASS pen needle boxes.

Please note that some products may not be available in all regions. Check with your local Terumo Medical Care Solutions - Pharmaceutical contact for details. All brand names are trademarks or registered trademarks owned by TERUMO CORPORATION, its affiliates, or unrelated third parties.



Terumo Europe N.V. +32 16 38 12 11
Terumo Corporation +81 3 6742 8500

Terumo Medical Corporation
Phone: +1 732 302 4900
Toll free (U.S.A.): +1 800 283 7866

Contact details are subject to change.
Please consult our website.

www.terumopharmaceuticalsolutions.com

Published by Terumo Europe 2023 - PM-06-441

Quality time for better care



PHARMACEUTICAL

Needle for pen-injectors
NANOPASS™ 32.5G



One size made to fit all

More quality of life for every patient

One size made to fit all

- The NANOPASS™ 32.5G range offers three different lengths – allowing you to choose the suitable pen needle for each patient type.
- Terumo has put its knowledge about reducing injection pain into the design of NANOPASS™ 32.5G. Users of long needles will also benefit from it.
- Color coding for easy identification of the preferred pen needle length.



NANOPASS™ 32.5G offers patients a more convenient insulin injection by applying dedicated Terumo technologies:

- 1 Asymmetric bevel¹**
- 2 Special thin pen needle size**
- 3 Terumo-thin wall technology**

For every type of patient the right needle!




Targeting on less pain², every day

1


Asymmetric bevel¹

The unique asymmetric bevel¹ of NANOPASS™ cuts the skin instead of puncturing it. This decreases the resistance while penetrating the skin – aims to reduce pain.²

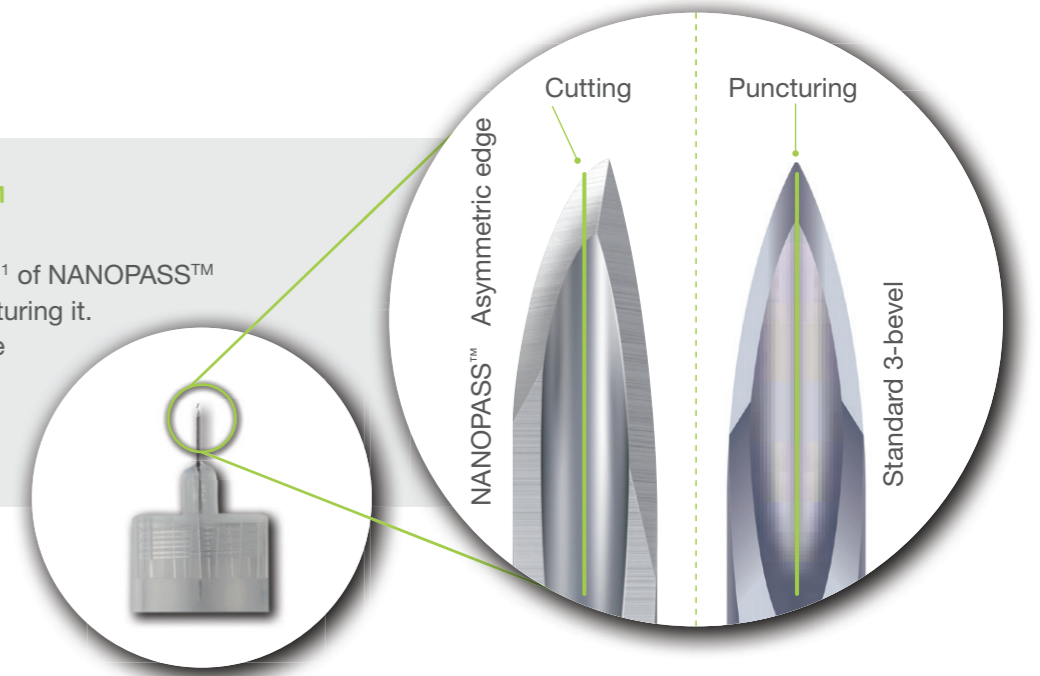

2

Special thin pen needle size

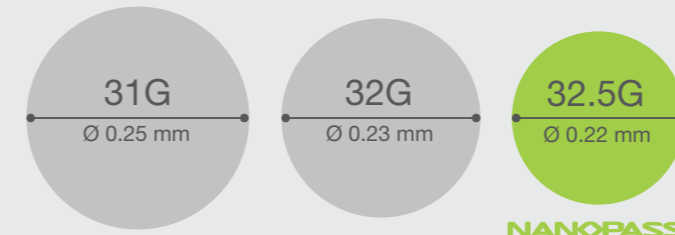

The special thin pen needle size of only 32.5G (0.22 mm) causes less skin trauma than larger needle sizes³.

3

Terumo-thin wall technology

The larger inner diameter of NANOPASS™ 32.5G optimizes the insulin flow for an easy insulin injection.

