



Quality time for better care

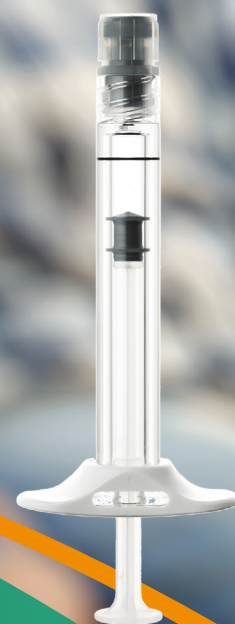

PHARMACEUTICAL

We've sharpened our focus on Ophthalmics

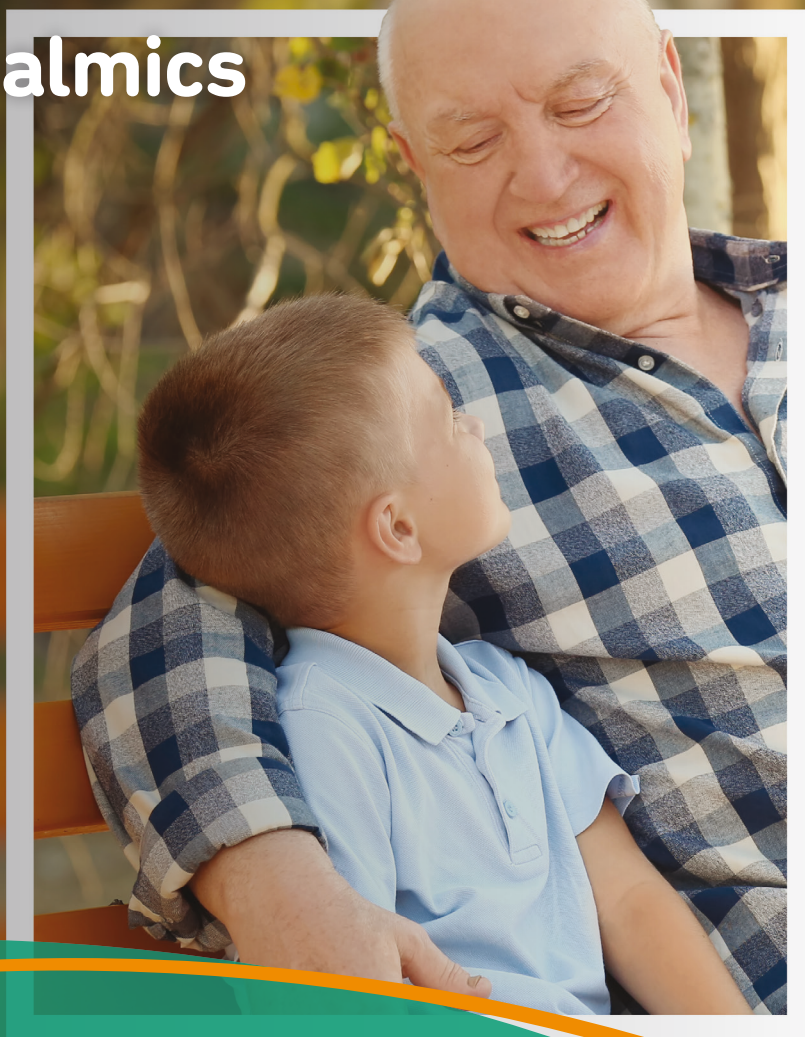

PLAJEX™ 0.5 ml Luer Lock silicone oil-free prefillable syringe

- Optimized for patient safety
- Ensuring dose accuracy
- Enhanced user experience

Integrated solution for challenging low dose applications

Designed to meet your needs and those of your patients, offering a best-in-class ready-to-fill polymer syringe solution for challenging low-dose applications, such as intravitreal injection drugs.

Optimized for patient safety

- Silicone oil-free syringe system enabled by proprietary i-coating™ technology
- Low subvisible particle level enables your drug product to fulfil the USP <789> requirements**

Ensuring dose accuracy

- Clear dose mark for dose confirmation
- Tight dimensional tolerances
- Easy-to-handle design

Enhanced user experience

- i-coating™ technology on stoppers for consistent and predictable break-loose and glide forces*
- Finger Grip with ergonomic design and backstop function
- Easy grip twist-off tip cap and integrated luer lock for secure needle connection



Specification	Description
Syringe barrel material	Cyclo Olefin Polymer (COP)
Syringe size	0.5 ml Reference to ISO 11040-6: 2019
Syringe barrel lubrication	None: Silicone oil-free
Dose mark	Ink-printed dose mark
Luer lock	Integrated Luer Lock tip EN ISO 80369-7
Flange shape	Cut flange Reference to ISO 11040-6: 2019
Twist-off tip cap • Rubber component • Rigid component	Chlorobutyl rubber (CIIR) Polypropylene
Plunger stopper	i-coating™ stopper Reference to ISO 11040-5: 2012 Chlorobutyl rubber (CIIR) with silicone resin coating
Packaging	Barrel: 10 x 10 Nest/Tub format (100 syringes per tub) Plunger stopper: Double Tyvek® bag or Rapid Transfer Port bag
Sterilization	Barrel: Autoclave (Steam) Plunger stopper: Autoclave (Steam)
Plunger rod	Polypropylene
Finger grip	Finger grip extension with backstop function, Acrylonitrile Butadiene Styrene (ABS) polymer

*Data on file

**Final drug product specifications are responsibility of pharmaceutical customers.

Images are for visualization purpose only.

Please note that some products may not be available in all regions, check with your local Terumo Medical Care Solutions - Pharmaceutical contact for details. All brand names are trademarks or registered trademarks owned by TERUMO CORPORATION, its affiliates, or unrelated third parties.

PLAJEX[†]
Ready-to-Fill Polymer Syringe

TERUMO

Terumo Europe N.V. **Terumo Corporation**
+32 16 38 12 11 +81 3 6742 8500

Terumo Medical Corporation
Phone: +1 732 302 4900
Toll free (U.S.A.): +1 800 283 7866

Contact details are subject to change.
Please consult our website.

www.terumopharmaceuticalsolutions.com